


Meditech Holter ECG monitors



Related products

- ❑ **CardioMera** and its Bluetooth-version, **CardioBlue24** can record continuous, full disclosure ECG up to 96 hours.

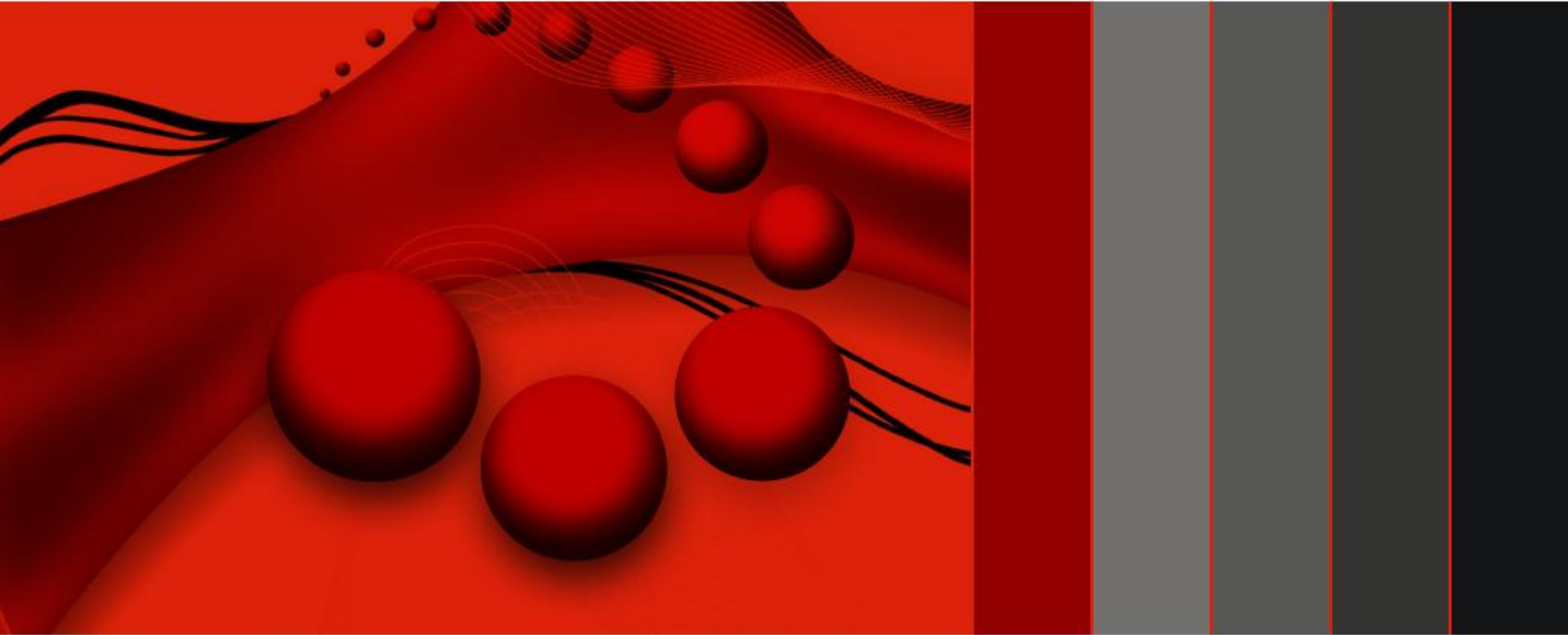

Indications for ambulatory ECG monitoring

- ❑ Arrhythmias
- ❑ Chest pain
- ❑ Unexpected syncope
- ❑ Shortness of breath
- ❑ Palpitations
- ❑ Evaluation of a pacemaker
- ❑ Regulation of anti-arrhythmic drugs
- ❑ Evaluation of patients after myocardial infarction
- ❑ Family history of heart disease

CardioMera/CardioBlue24








hardware features

- ❑ Dimensions: 66*59x17 mm
- ❑ Weight: 50 g
- ❑ Sampling frequency: 1200 or 600 Hz
- ❑ ECG storage: 150, 300 or 600 Hz on a standard memory card
- ❑ Storage capacity: up to 96 hours of ECG
- ❑ ECG recording: 1, 2, 3, 5 channels (7 leads)
- ❑ Data transfer: USB optical cable (CardioMera) or Bluetooth (CardioBlue24)
- ❑ Power consumption: 1 AAA alkaline or rechargeable battery
- ❑ Operator control: 1 pushbutton
- ❑ Warranty: 2 years for the main recorder unit

CardioVisions software features

- ❑ Integrated patient database
- ❑ Personal or network use
- ❑ Automatic backup function for easy restore of data
- ❑ Online or offline programming: the unit can be programmed either by connecting the device to the PC or by inserting a pre-configured flash card into the device
- ❑ Real-time ECG monitoring for better electrode placement
- ❑ Full disclosure Holter ECG annotation
- ❑ Manual annotation editing for more accurate reports
- ❑ Time and frequency domain HRV
- ❑ Pacemaker function
- ❑ QT dispersion
- ❑ ST and HRV statistics
- ❑ Graphical, tabular and textual display of data
- ❑ E-mail function for recordings and reports
- ❑ Automatic report generator for faster evaluation
- ❑ Editable and printable reports

CardioVisions software outlook

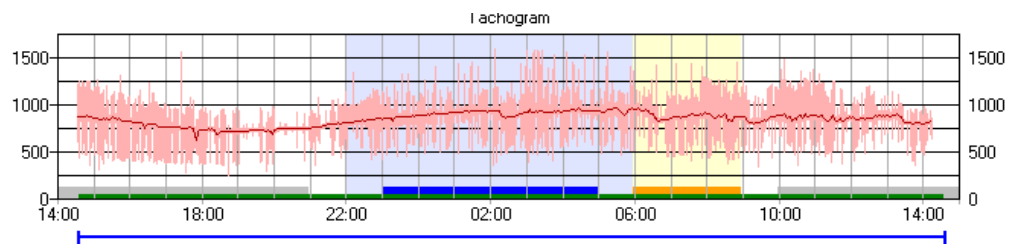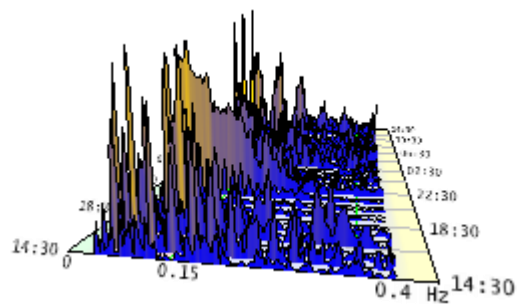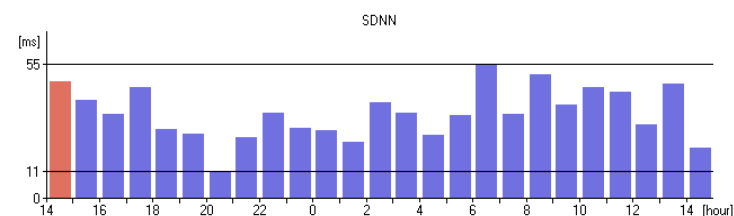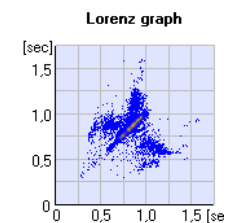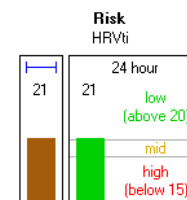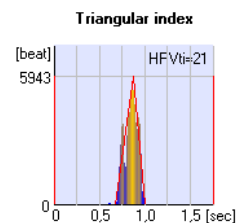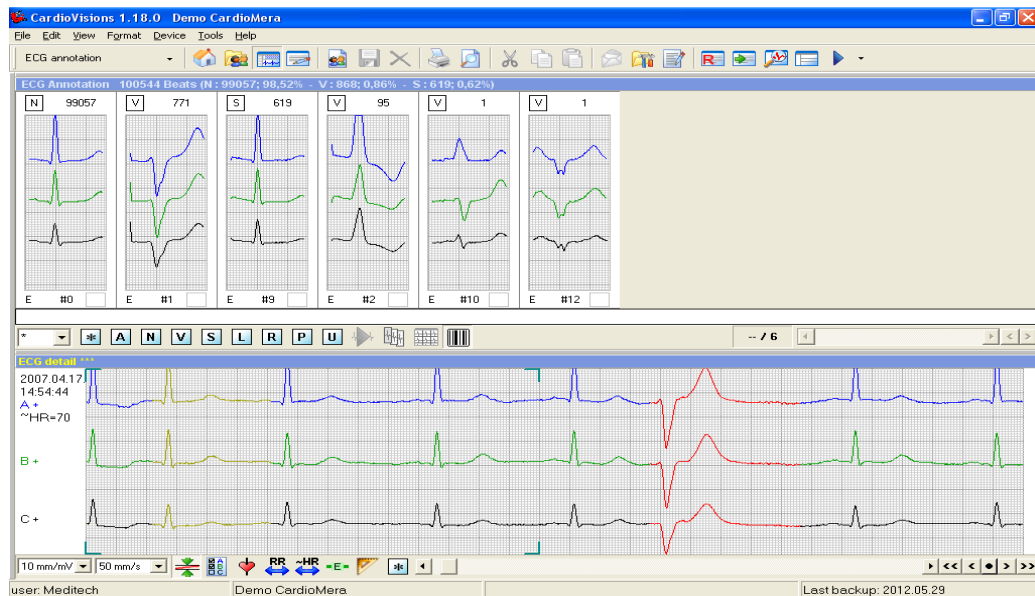

Set content

- 1* ECG recorder unit
- 1* memory card
- 1* USB cable (CardioMera) or Bluetooth USB dongle (CardioBlue24)
- 7*leads of patient ECG cable
- 1*package of disposable ECG electrodes
- 1* AAA accumulator
- 1* CD with the latest software and manual
- 1*user documentation (quality certificate)

Support

- ❑ Guide for software installation and registration
- ❑ Video guide for software installation and registration
- ❑ Guide for programming the ECG Holter monitors and evaluating recorded ECG

Contact details

For more information please contact:

Meditech Ltd.

1184 Budapest, Mikszáth Kálmán utca 24., Hungary

Phone: +36 1 280 8232, +36 1 280 8233

Fax: +36 1 282 9388

E-mail: meditech@meditech.eu

Web: www.meditech.eu

Thank you for your attention!